

RABGGTA Polyclonal Antibody
Purified Rabbit Polyclonal Antibody (Pab)
Catalog # AP57862**Specification**

RABGGTA Polyclonal Antibody - Product Information

Application	IHC-P, IHC-F, IF, ICC
Primary Accession	Q92696
Reactivity	Rat, Pig
Host	Rabbit
Clonality	Polyclonal
Calculated MW	65072

RABGGTA Polyclonal Antibody - Additional Information**Gene ID** 5875**Other Names**

Geranylgeranyl transferase type-2 subunit alpha, 2.5.1.60, Geranylgeranyl transferase type II subunit alpha, Rab geranyl-geranyltransferase subunit alpha, Rab GG transferase alpha, Rab GGTase alpha, Rab geranylgeranyltransferase subunit alpha, RABGGTA

Dilution

IHC-P ~ ~ N/A
IHC-F ~ ~ N/A
IF ~ ~ 1:50 ~ 200
ICC ~ ~ N/A

Format

0.01M TBS(pH7.4), 0.09% (W/V) sodium azide and 50% Glyce

Storage

Store at -20 °C for one year. Avoid repeated freeze/thaw cycles. When reconstituted in sterile pH 7.4 0.01M PBS or diluent of antibody the antibody is stable for at least two weeks at 2-4 °C.

RABGGTA Polyclonal Antibody - Protein Information**Name** RABGGTA**Function**


Catalyzes the transfer of a geranylgeranyl moiety from geranylgeranyl diphosphate to both cysteines of Rab proteins with the C-terminal sequence -XXCC, -XCXC and -CCXX, such as RAB1A, RAB3A, RAB5A and RAB7A.

RABGGTA Polyclonal Antibody - Protocols

Provided below are standard protocols that you may find useful for product applications.

- [Western Blot](#)
- [Blocking Peptides](#)
- [Dot Blot](#)
- [Immunohistochemistry](#)
- [Immunofluorescence](#)
- [Immunoprecipitation](#)
- [Flow Cytometry](#)
- [Cell Culture](#)

RABGGTA Polyclonal Antibody - Images

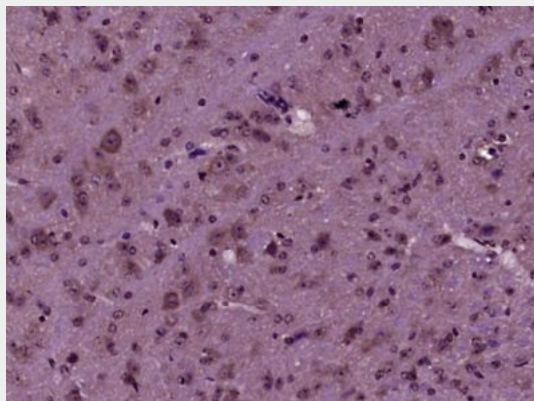

Paraformaldehyde-fixed, paraffin embedded (Mouse brain); Antigen retrieval by boiling in sodium citrate buffer (pH6.0) for 15min; Block endogenous peroxidase by 3% hydrogen peroxide for 20 minutes; Blocking buffer (normal goat serum) at 37°C for 30min; Antibody incubation with (RABGGTA) Polyclonal Antibody, Unconjugated (bs-21088R) at 1:400 overnight at 4°C, followed by operating according to SP Kit(Rabbit) (sp-0023) instructions and DAB staining.